

PCR Plates

96 and 384-well, standard, half and full skirted




96-well Non-Skirted 200 µl



96-well Half-Skirted 200 µl



96-well Full-Skirted 200 µl

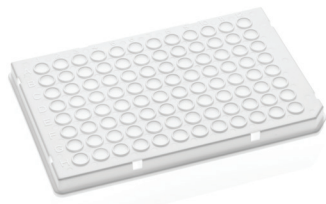

96-well for LightCycler™ 100 µl



Full Skirted 384-well 50 µl


PRODUCT HIGHLIGHTS

Developed for use with the industry's most popular thermocyclers, Biotix PCR plates are manufactured and tested to ensure robust and reproducible amplification under the most exacting protocols. Made from medical-grade virgin polypropylene, these plates feature ultra-thin wall wells that provide excellent thermal transfer with PCR heat blocks. Each plate offers a molded-in alphanumeric grid to help with sample


identification, and can withstand centrifugation up to 4,000 RCF when supported by the appropriate plate rotor adaptor. Slightly raised well rims facilitate secure sealing with films or strip caps. As with all Biotix products, PCR plates are free from heavy metals, mold releases or other extractables that could interfere with highly sensitive protocols. They are BioReady™ Certified for purity.

FEATURES	BENEFITS
Ultra-clear, medical-grade virgin polypropylene	Biologically inert, offers excellent chemical and temperature resistance with superior clarity for sample visualization
Ultra-thin walled construction	Ensures excellent thermal transfer with PCR heat blocks and uniform amplification
Available in standard, half-skirt or full skirt	Broad range of compatibility with leading thermocycler brands; half and full-skirted plates enable labeling on sides of plate
BioReady™ Certified	PCR plates are certified free of RNase, DNase and endotoxins (pyrogens)


Product Dimensions

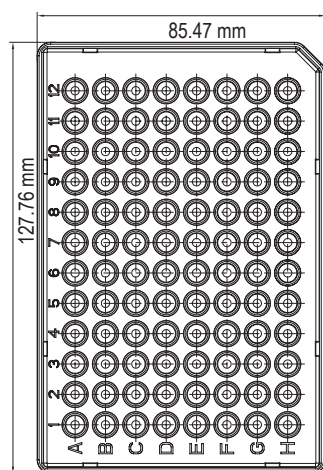


96-well Half-Skirted 200 µl

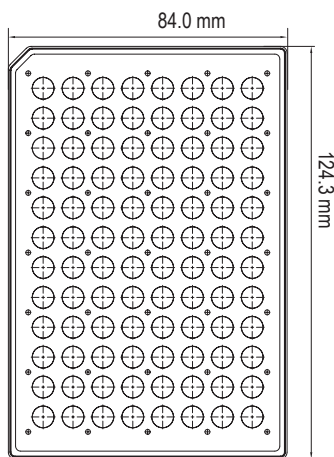


96-well Non-Skirted 200 µl

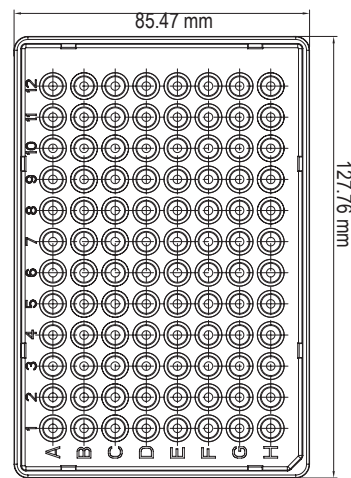

96-well Full-Skirted 200 µl

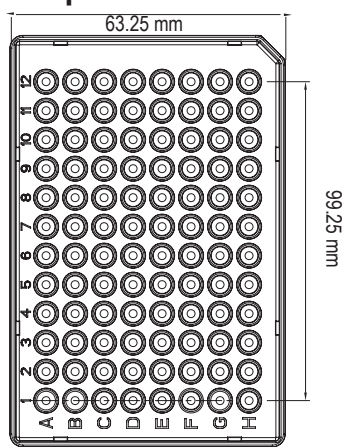

96-well for LightCycler™ 100 µl

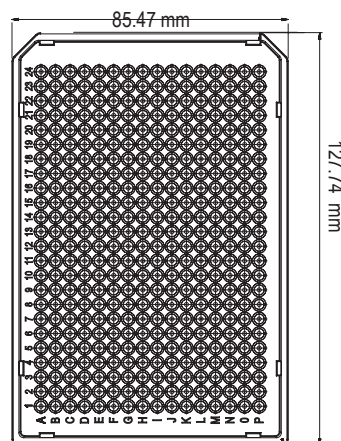

384-well Full Skirted 50 µl

Ordering Information

Biotix Part No.	Description	Volume (Max. capacity/well)	Packaging
P-0200-9C	PCR Plate, 96-well, Non-Skirted, Clear	200µl	10 plates/pack 5 packs/case
P-0200-9CHS	PCR Plate, 96-well, Half-Skirted, Clear	200µl	10 plates/pack 5 packs/case
P-0200-9CFS	PCR Plate, 96-well, Full-Skirted, Clear	200µl	10 plates/pack 5 packs/case
P-0100-9CHS	PCR Plate, 96-well, for LightCycler™, White	100µl	10 plates/pack 5 packs/case
P-0050-3CFS	PCR Plate, 384-well, Full-Skirted, Clear	50µl	10 plates/pack 5 packs/case

©2013, Biotix, Inc. All rights reserved

Biotix, Fluid Innovation, Delta Filter, X-Resin, Blade, FlexFit, Statix, Leave Nothing Behind, BioReady, Tip Eject, and XTreme Quality are registered trademarks or trademarks of Biotix, Inc. All other brands and names contained herein are the property of their respective owners.